

Key Features

- Spiked drip chamber with integral air inlet**
 Suitable for use with non-collapsing fluid containers (eg glass) without the expense or needle-stick risk of using a needle mounted air inlet
- Integral 15mm disc filter in drip chamber**
 Prevents accidental administration of particulate matter in infused fluid
- Roller clamp distal to drip chamber**
 Accurate control of fluid infusion rate at the tip of your thumb
- Built in needle free valve (GNF sets only)**
 Permits administration of intermittent bolus dose of drugs without the risks associated with administration via a needle or a 3-way stop tap
- Luer thread on female receiver of needle free valve**
 Luer lock syringes can be screw locked onto needle free valve
- Very low DEHP plastics**
 Administration set tubing far exceeds European safety standards for DEHP containing plastics (European standard maximum – 0.1%, AMS G-Sets – 0.003%)


GS1600-Set

Product code – GS1600

Gravity infusion set

L – 160 cm, prime vol – 12 ml

Units / box - 200

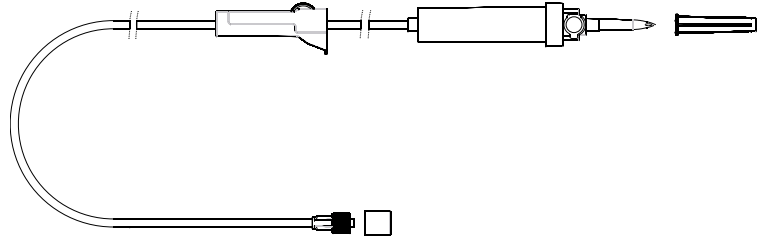


GNF2000-Set

Product code – GNF2000

Gravity infusion set with integral needle free valve

L – 200 cm, prime vol – 15 ml

Units / box - 200

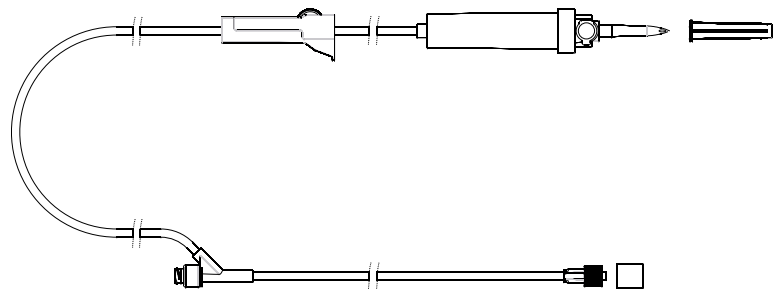


GNF3500-Set

Product code – GNF3500

Gravity infusion set with integral needle free valve

L – 350 cm, prime vol – 25 ml

Units / box - 200

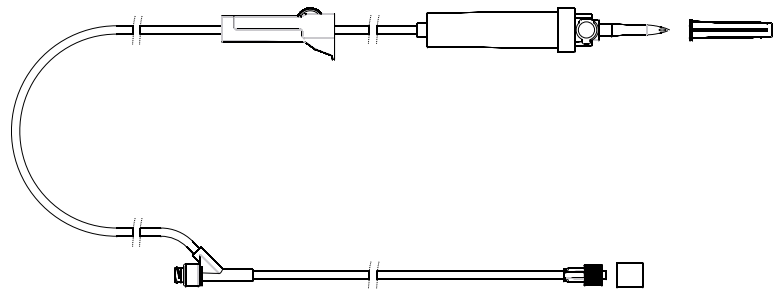

Accessories

STNF100

Product code – STNF100

3-Way stop tap with needle free valve and extension

L – 100 cm, prime vol – 7.1 ml

Units / box 100

